



VACCINE PREVENTABLE DISEASES REFERENCE CENTERS TESTING RESOURCES

The Vaccine Preventable Diseases Reference Centers (VPD RCs) network is composed of four public health laboratories that provide diagnostic testing and typing for seven vaccine preventable diseases (Measles, Mumps, Rubella, Varicella-zoster, *B. pertussis*, *N. meningitidis* and *H. influenzae*) as well as Middle East respiratory syndrome (MERS-CoV). Since July 2023, VPD RCs also support [enterovirus surveillance](#) through screening and typing.

Testing is performed using standardized methods developed by CDC and is available to public health departments free of charge. Submitting sites are assigned to one or two of the VPD RCs depending on the services requested.

Figure 1. VPD RC Assignments
Visit aphl.org/VPD for the most up-to-date assignments

VPD RC Laboratories


Viral Pathogen Testing

- California Department of Public Health Laboratory (Contact: VRDL.submittal@cdph.ca.gov)
- New York State Department of Health: Wadsworth Center (Contact: VPDNYS@health.ny.gov)
- Minnesota Public Health Laboratory Division (Contact: health.vpd.submissions@state.mn.us)
- Wisconsin State Laboratory of Hygiene (Contact: cddvpdproject@slh.wisc.edu)

Bacterial Pathogen Testing

- Minnesota Public Health Laboratory Division
- Wisconsin State Laboratory of Hygiene

Testing for *Streptococcus pneumoniae* is available through the [Antimicrobial Resistance Network \(AR Lab Network\)](#).

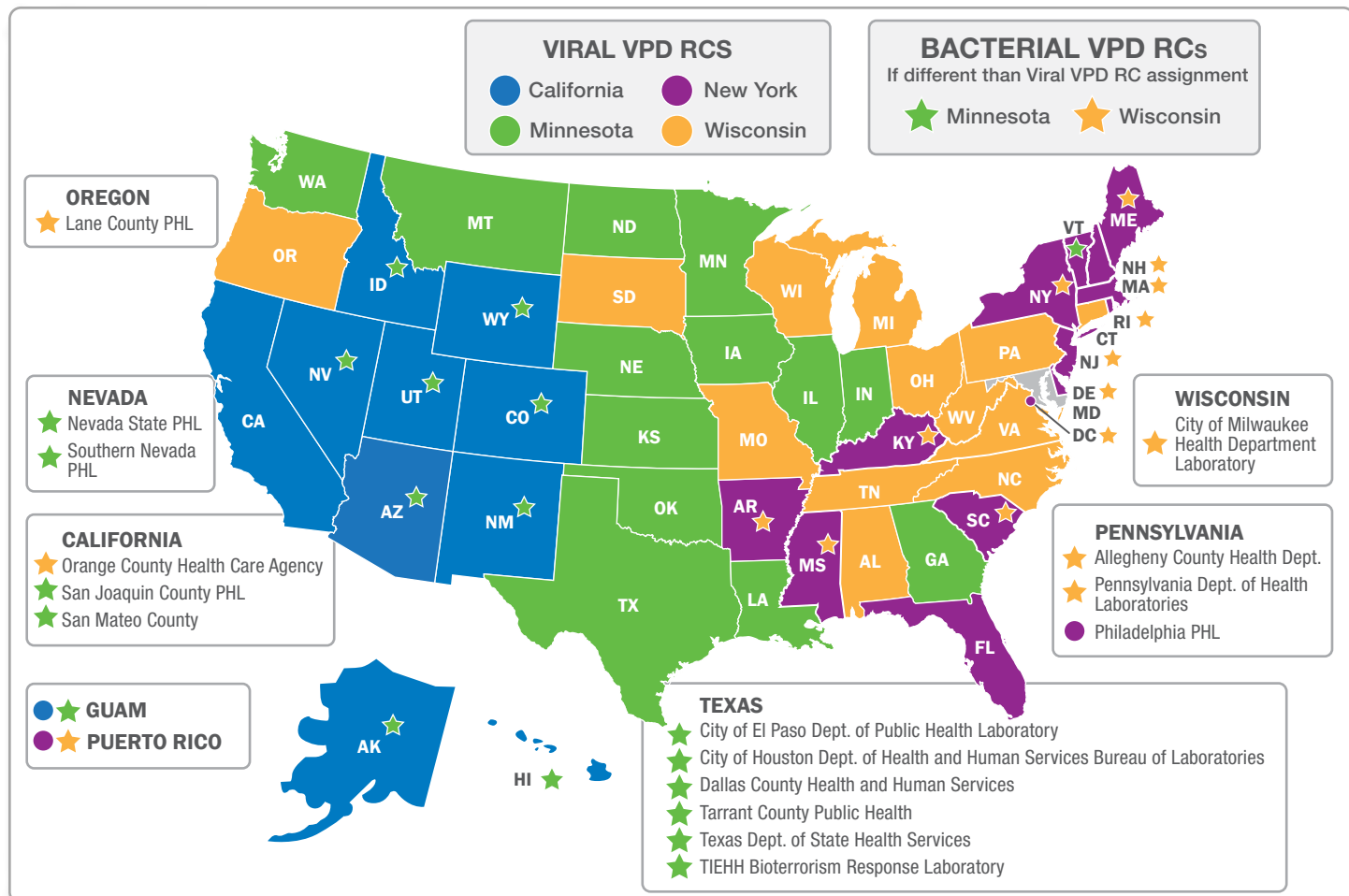

Testing Available at VPD RCs

Bacterial Pathogen Testing

VPD RC laboratories that provide bacterial pathogen testing perform: PCR for detection, serotyping/serogrouping and culture.

Pathogen	PCR	Serotyping/ Grouping	Maximum Turnaround Time (# of Business Days)		Location(s)/ Network
			PCR	Serotyping/ Grouping	
<i>B. pertussis</i> *	Yes	N/A	2	N/A	Minnesota Only
<i>H. influenzae</i>	Yes	Yes	2	5	Both Bacterial RCs
<i>N. meningitidis</i>	Yes	Yes	2	5	Both Bacterial RCs
<i>S. pneumoniae</i>	Yes	Yes	2	5	AR Lab Network†

Viral Pathogen Testing

VPD RC laboratories that provide viral pathogen testing perform PCR and genotyping, which is conducted on all PCR-positive specimens unless otherwise indicated as part of larger outbreak.

Pathogen	PCR	Genotyping	Maximum Turnaround Time (# of Business Days)		Location(s)
			PCR	Genotyping	
Measles virus‡	Yes	Yes	2	10	All Viral RCs
Mumps virus	Yes	Yes	2	10	All Viral RCs
Rubella virus	Yes	Yes	2	10	All Viral RCs
Varicella-zoster virus	Yes	Yes	2	10	All Viral RCs
Enterovirus	Yes	Yes	Surveillance testing only, results not reported		All Viral RCs
MERS-CoV	Yes	N/A	2	N/A	California & New York Only

* Culture of *B. pertussis* is also available with a turnaround time of 10 business days.

† Contact Minnesota (arlnmn@state.mn.us) or Wisconsin (wIARLN@slh.wisc.edu) for current AR Lab Network submission instructions.

‡ Measles virus vaccine assay (MeVA) has a maximum turnaround time of three business days.

Reporting Results from VPD RCs

The VPD RCs will receive specimens and perform the appropriate tests,[§] with anticipated turnaround times as listed in the testing tables on Page 2. Results will be reported, with patient identifiers, by the RC to the submitting laboratory through secure web portal, encrypted email, secure fax or over the phone. Results for all specimens tested will be reported to CDC by the RC via electronic HL7 messaging or a secure FTP site; however, these data will not replace state notifications to CDC through the Nationally Notifiable Diseases Surveillance System, which should be submitted per the usual protocol.

VPD RC Specimen and Shipping Recommendations

It is critical that specimens are properly packaged and shipped to the VPD RCs to ensure optimal test results and avoid delays in testing or rejection of specimens. Please contact your assigned VPD RC (**Figure 1**) with questions about specimen submission and to obtain specimen requisition information/forms. Submitting agencies should phone or email to notify the RC about intention to submit eligible specimens. Shipping costs are covered by VPD funding; please contact your VPD RC or APHL to obtain that information. See below for specimen submission requirements for available VPD RC viral and bacterial pathogen testing.

Bacterial Pathogen Recommendations

Note: Swabbed specimens must be collected with **synthetic swabs**; cotton swabs are **not acceptable**.

Method	Specimen Type	Minimum Volume	Specimen Storage	Shipping Recommendations
<i>B. pertussis</i> Real time PCR	NP swab or isolate	600 µL	<ul style="list-style-type: none"> • Swabs: Refrigerate at 4 °C as soon as possible. • Isolates: Refrigerate in Regan-Lowe transport medium or freeze on cryobeads. 	<ul style="list-style-type: none"> • Refrigerated specimens: Ship on cold packs. • Frozen specimens: Ship on dry ice.
<i>B. pertussis</i> Culture	NP aspiration, wash or swab	250 µL	<p>Samples must be received by the RC within 72 hours of collection.</p> <ul style="list-style-type: none"> • If plated within 24 hours of collection: Use Amies transport medium with charcoal. • For delays of 24–72 hours: Use Regan-Lowe transport medium (available from MDH-PHL). 	<p>All media specimens:</p> <ul style="list-style-type: none"> • Ship overnight. • Ambient temperature.
<i>N. meningitidis</i> PCR, Serogrouping	CSF or isolate	500 µL	<ul style="list-style-type: none"> • Primary specimen: Freeze. • Isolates: Transport on chocolate agar slants at ambient temperature. 	<ul style="list-style-type: none"> • Frozen specimens: Ship on dry ice. • Isolates: Ship at ambient temperature.
<i>H. influenzae</i>	CSF or isolate	500 µL		
<i>S. pneumoniae</i> Real time PCR, serotyping	CSF or isolate	250 µL		

§ Additional testing may be performed on submitted specimens.

Viral Pathogen Recommendations

Note: Swabbed specimens must be collected with **synthetic swabs** in **viral transport media (VTM/UTM) only**; cotton swabs and molecular transport media (MTM) are **not acceptable**.

Method	Specimen Type	Minimum Volume	Shipping Recommendations
			All specimens types: Ship overnight, Monday–Thursday
Measles virus PCR	<ul style="list-style-type: none"> Throat swab in VTM/UTM Nasopharyngeal (NP) swab in VTM/UTM Combined throat/NP swab in VTM/UTM Urine Cerebrospinal fluid (CSF) (MN only, in cases of encephalitis) 	250 µL	<p>Storage: Place swabs in 2 mL standard VTM/UTM.</p> <p>Shipping:</p> <ul style="list-style-type: none"> If shipping within 24 hours: Store at 4°C If shipping is delayed: Store at -70°C, ship on dry ice
	<ul style="list-style-type: none"> Buccal swab in VTM/UTM NP swab in VTM/UTM (not WI) Throat swab in VTM/UTM (not WI) Urine (MN only) CSF (MN only, in cases of encephalitis) 	250 µL	
Rubella virus PCR	<ul style="list-style-type: none"> Throat swab in VTM/UTM NP swab in VTM/UTM Urine (MN & WI only) CSF (MN only, in cases of encephalitis) 	250 µL	
Enterovirus	<ul style="list-style-type: none"> NP swab in VTM/UTM CSF 	300 µL	
	Stool	Solid: 1 g Liquid: 300 µL	
Varicella-zoster virus PCR	Scab	N/A	
	Skin lesion swab in VTM/UTM	250 µL	Ship at 4°C or frozen on dry ice.
	CSF	250 µL	Freeze, ship on dry ice.
MERS-CoV PCR	<ul style="list-style-type: none"> Lower respiratory specimens Upper respiratory swabs in VTM/UTM Serum 	250 µL	<ul style="list-style-type: none"> If collected within three days of receipt at RC: Refrigerate, ship on frozen gel packs. If collected over three days prior to receipt at RC: Freeze, ship on dry ice.
	Stool (CA only)	Solid: 1 g Liquid: 300 µL	